


TOPAZ[®]

Desktop Video Magnifier

Features at a Glance

- Clear, vivid images
- Magnifies up to 77x
- Freeze Frame, Focus Lock, Find feature, and Position Locator
- Directed LED lighting for a bright work surface without glare
- 8.25 inches of workspace under the camera
- 30 high-contrast text color modes plus true color and enhanced full color
- Smooth XY reading table
- Monitor raises, lowers, tilts, and swivels
- Unique groove on XY table holds bottles and cans still for viewing

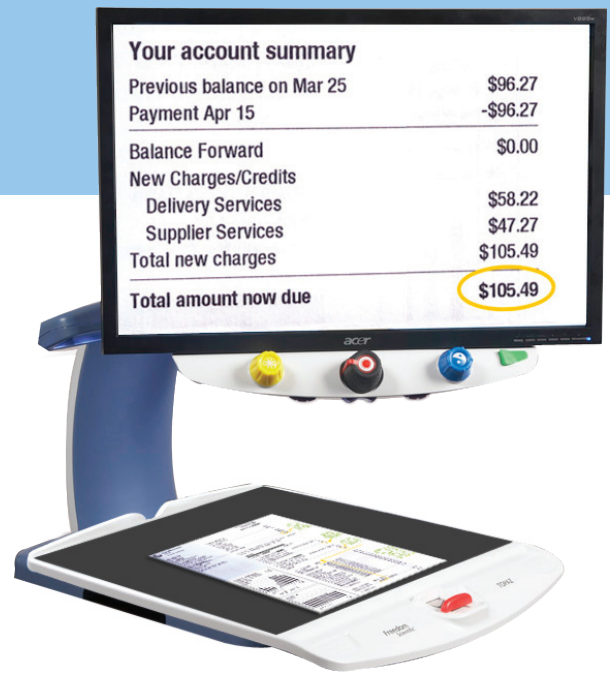

Versatile, High-Powered Magnification for Home, School or Office

If low vision keeps you from reading and enjoying your favorite hobbies and crafts, the TOPAZ desktop video magnifier can help. With magnification all the way up to 77x, right away you'll be reading your mail, clearly seeing prescription instructions, and more. With a large reading table and 8.25 inches between the camera and the table, there's also plenty of room for hobbies and crafts.

The attractive TOPAZ fits in any environment – home, school, or office, and ships with your choice of monitor option: a standard 17-inch monitor or 20-inch, 22-inch, or 24-inch widescreen. You can also choose our camera only model to use the TOPAZ with your own monitor instead.

Features Others Don't Have

Focus Lock keeps the object or document you are magnifying crisp and clear when you're writing under the camera. The **Find Feature** zooms out so you can quickly reposition your document and zoom back in to your working magnification. **Freeze Frame** is great for referencing items from a list or for magnifying information when you need your hands free for other activities.



Learn more online at www.FreedomScientific.com

TOPAZ[®]

Desktop Video Magnifier

Magnification:

17-inch monitor - 1.9x - 49x

20-inch monitor - 2.4x - 64x

22-inch monitor - 2.6x - 68x

24-inch monitor - 2.9x - 77x

Color Modes:

Full Color, Black on White, White on Black, Yellow on Blue, Yellow on Black, plus Enhanced Full Color and 27 additional, configurable color modes

Working Space Under Camera:

8.25 in. (21 cm.)

Dimensions:

25 in. x 17.2 in. x 17.3 in.

(63.5 cm. x 43.7 cm. x 43.9 cm.)

Total Unit Weight:

Camera Only: 32.8 lbs. (14.1 kg.)

17-inch monitor - 39.5 lbs. (17.9 kg.)

20-inch monitor - 41 lbs. (18.6 kg.)

22-inch monitor - 42.1 lbs. (19.1 kg.)

24-inch monitor - 43.7 lbs. (19.8 kg.)

Controls:

On / Off toggle switch

Color mode select dial

Magnification 16-position dial

Freeze Frame push button

Brightness dial

Find push button

Focus Lock push button

Power Supply:

Single wall plug, 50-60 Hz , 100-240V AC



Freeze Frame feature keeps your hands free for other activities

Designed & Assembled
in the USA 



**Freedom
Scientific[®]**
everything in sight.[™]

Call to Learn More
800-444-4443

www.FreedomScientific.com


*Freedom Scientific is the world's leading
manufacturer of assistive technology products for
those who are low vision or blind.*

Visit us online at **www.FreedomScientific.com** to learn more about these and other innovative products from Freedom Scientific that make living with low vision easier.



RUBY[®]

Pocket-Sized Handheld Video
Magnifier for maximum portability

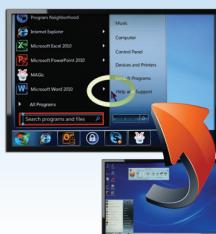

MAGic[®]

HD Screen Magnification Software
with built-in screen reading features



ONYX[®] Deskset XL

Three-way Portable Video Magnifier
for easy magnification near or far