


# TOPAZ<sup>®</sup> EZ HD Control Panel Reference Card

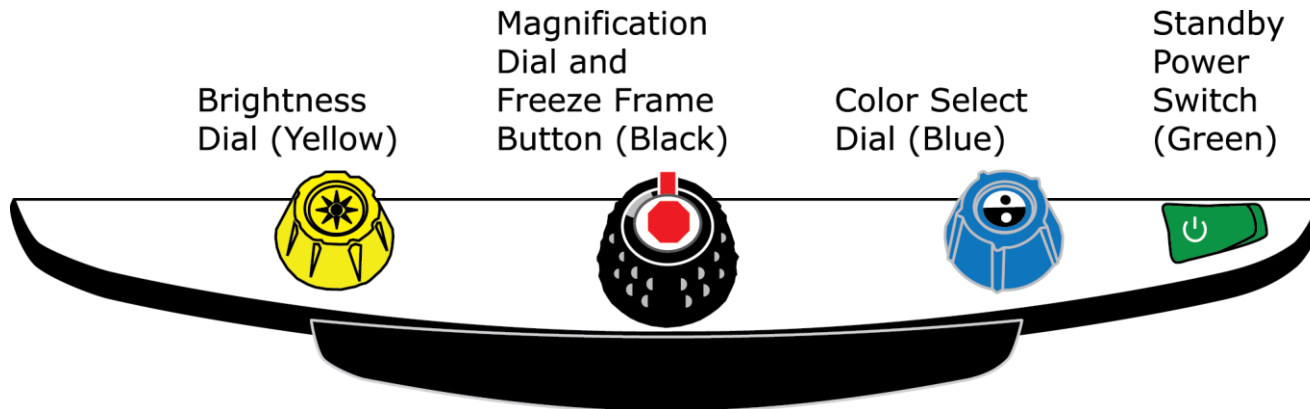

Dial / Button / Switch	Functionality
Brightness Dial (Yellow)	<ul style="list-style-type: none"><li>• Turn to increase or decrease screen brightness</li></ul>
Magnification Dial / Freeze Frame Button (Black)	<ul style="list-style-type: none"><li>• Turn to increase or decrease magnification</li><li>• Press center to switch Freeze Frame on or off</li></ul>
Color Select Dial (Blue)	<ul style="list-style-type: none"><li>• Turn to choose a color mode that best suits your vision needs</li></ul>
Standby Power Switch (Green)	<ul style="list-style-type: none"><li>• Press to the right to place unit in standby mode</li><li>• Press to the left to turn on unit and exit standby mode</li></ul>